HORICO DENTAL

100
1918 - 2018

HOPF, RINGLEB & CO. GMBH & CIE.

Alt. Nummernsystem - ALT. NUMBERING SYSTEM
SYSTEME DE NUMÉRATION ALT. - ALT. SISTEMA DE NUMERACIÓN

HORICO Bestell-Nr.
HORICO ORDER NO.

Diamantschleifer
DIAMOND INSTRUMENTS

Form	HORICO	Seite - PAGE
SHAPE	Form - SHAPE	FG HP/PM
368	257	29 55
379	277	30 55
389	494	16
390	272	29 55
392	465	32
515	515	32
801	001	15 51
801L	697	15
801M	001L	15
802	002	16 51
805	010/014	17 51
806	019	17 52
807/822	225	17 51
811	038	18
812	023	52
814	023	52
818	041	18 52
820	042	52
824A	055	18
825	304	32 56
830	233	26
830A	239	27 54
830L	234	26
830R	237	27
830RL	238	27
833	466	32
834	834	33
835	108/109	18/19 52
835A	108	18
835KR	156	22
836	110	19
836KR	157	22
837	111	19 52
837KR	158	23
837L	112	19 52
837LKR	158L	23
838	138/139	21
838A	138	21
838M	138L	21
839	150	22
839K	150K	22
842	113	19 52
842R	143	46 53
845	170	24
845KR	544	24
846	171	24
846KR	545	24
847	172	24 53
847KR	546	24
848	173	25 54
848L	174	25
848LE	174	34
850	199	25 54
850L	200	26
851	218	26
852	164	23 53
855	197	25
855LE	200	35
856	198	25 54
857	220	26
857LE	188	34
858	165	23 53
859	166	23 53
859L	167	23 53
860	247	28 54
861	248	28 54
862	273	29
863	274	29 55
863LE	256	35
864	251	28 55
875	535	32
877	288	30
877K	297	31
878	289	30
878K	298	31
878P	255A	28
879	290	30
879K	299	31
879L	291	31 55
879P	299A	31
880	140	21
881	141	21 53
881P	141A	21
882	142	22 53
883	539	32 56
884	129	20
885	130	20

Diamantschleifer
DIAMOND INSTRUMENTS

Form	HORICO	Seite - PAGE
SHAPE	Form - SHAPE	FG HP/PM
886	131	20
888	249	28 54
888L	250	28 55
888LE	250	28
889	540	32 56
892	292	55
899	033	17
899E	033	17
909	068	18 52
955	699	32
956	164	23 53

Hartmetall-Fräser
TUNGSTEN CARBIDE CUTTERS

Form	Verzahnung	HORICO	Seite - PAGE
SHAPE	TOOTHING	SHAPE TOOTHING	PAGE - PAGINA
C129	E	141 190	81
C129	E	140 190	81
C129	EF	141 140	85
C129	FSQ	141 134	88
C129	NE	141 NE	89
C129	NEF	141 NEF	89
C136	E	197 190	82
C136	EF	197 140	85
C136	EVF	197 110	-
C137	E	225 190	82
C137	EF	225 140	85
C137KR	E	235 190	-
C137KR	EF	235 140	-
C138	E	198 190	82
C138	EF	198 140	85
C138	EVF	198 110	90
C138	FSQ	198 134	88
C139	DF	289 141	-
C139	E	289 190	83
C139	EF	289 140	86
C139	EUF	289 110	90
C139	FSQ	289 134	88
C139	FSQL	289 132	90
C139	NE	289 NE	89
C139	NEF	289 NEF	89
C219		468 211	81
C250	AQ	277 137	87
C250	E	277 190	83
C250	EF	277 140	86
C250	FSQ	277 134	88
C251		274 175	84
C251	AQ	274 137	87
C251	AQL	274 139	90
C251	E	274 190	83
C251	EAL	274 192	90
C251	EF	274 140	86
C251	FSQ	274 134	88
C251	FSQA	274 136	87
C251	GE	274 220	80
C251	GSQ	274 212	81
C251	KX	274 KX	89
C251	SGE	274 223	80
C251	SGEA	274 225	80
C257	E	187 190	81
C257	EF	187 140	85
C257	Q	161 176	83
C257R		200 175	84
C257R	E	200 190	82
C257R	EF	200 140	85
C257R	FSQ	200 134	88
C261		194 175	84
C261	AQ	194 137	87
C261	E	194 190	82
C261	EF	194 140	85
C261	EVF	194 110	90
C261	FSQ	194 134	88
C261	FSQA	194 136	87
C261	FSQL	194 132	90
C261	GSQ	194 212	81
C351	E	263 190	83
C283	EF	289 140	86
C283	FSQ	289 134	88
C295	E	292 190	83
C295	EF	292 140	86
C295	FSQ	292 134	88
C296	E	110 190	81

Hartmetall-Fräser
TUNGSTEN CARBIDE CUTTERS

Form	Verzahnung	HORICO	Seite - PAGE
SHAPE	TOOTHING	SHAPE TOOTHING	PAGE - PAGINA
C296R	E	140 190	81
C296R	EF	140 140	84
C351	AQ	263 137	87
C351	E	198 190	82
C351	E	263 190	83
C351	EF	198 140	85
C351	EF	263 140	86
C351	FSQ	198 134	88
C356R	E	200 190	91
C364	E	110 190	81
C364R	E	144 190	81
C364R	EF	144 140	85
C364R	EVF	144 110	81
C42	EF	010 140	84
C71		001 175	84
C71	E	001 190	81
C71	EF	001 140	84
C71	FSQ	001 134	88
C72	AQ	144 137	87
C72	E	144 190	81
C72	EF	144 140	85
C72	FSQ	144 134	88
C72	GE	144 220	80
C73	E	277 190	83
C73	EF	277 140	86
C73	EUF	277 110	90
C73	FSQ	277 134	88
C73	NE	277 NE	89
C73	NEF	277 NEF	89
C76	E	254 190	82
C77	E	237 190	82
C77	E	239 190	82
C77	EF	237 140	86
C77	EF	239 140	86
C77	EVF	237 110	90
C77	FSQ	237 134	88
C78	E	254 190	82
C78	E	257 190	83
C78	EF	257 140	86
C79		194 175	84
C79	AQ	194 137	87
C79	E	194 190	82
C79	EF	194 140	85
C79	EF	194 138	86
C79	EVF	194 110	90
C79	FSQ	194 134	88
C79	FSQA	194 136	87
C79	FSQL	194 132	90
C79	GE	194 220	80
C79	GSQ	194 212	81
C79	KX	194 KX	89
C79	NE	194 NE	89
C79	NEF	194 NEF	89
C79	SGE	194 223	80
C79	SGEA	194 225	80
C79	SGEL	194 224	90
C85	EFAL	272 142	90
C85	FSQ	272 134	88
		113 137	78
		113 190	78
		113 212	78

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HORICO® offers you...

... top quality

All of our medical products are certified. And we take this even a step still further: Every single item that leaves any of our plants has been subjected to a detailed visual inspection under the microscope, and to other tests. However, our great manufacturing depth is the basis for our supreme quality. From the raw materials to the completed product, almost every step of production is carried out in our own facilities in Germany. All over the world, our brand HORICO® is recognized as a sign of superior quality and reliability.

... Innovation

From the start, we have been constantly developing new products to find the optimal solutions for our customers. Some solutions are still in use, such as our Steelcarbo® strips, which were the first product of our company. Other highlights of our research and development:

1918 - the first separating strips not made of paper (sandpaper) but of steel.

1930 - start of production of diamond burs, which today still forms our core business, thanks to their outstanding features.

1956 - Diaflex® Disc: Invention of the world's first flexible and stainless diamond disc with a thickness of .20 or .30 mm

1963 - Diatrepano: Uniquely shaped diamond tool for the processing of cavities and pulp chamber opening.

1974 - Superdiaflex® disc, even more flexible with a thickness as low as .10 or .15 mm

1984 - Superapid®: Development of the world's first slightly flexible sintered diamond disc, which will not break as easily as former discs.

1986 - Expansion of HORICO's program to include carbide products. By now, in the sector of carbide cutters we are also among the global leaders, offering all popular kinds of toothings.

1998 - Inclusion of a ball with a diameter of no more than .50 mm for ultra fine work into the program.

2008 - New surface treatment for carbide cutters. The new „black“ carbide cutters are highly successful because of their combination of excellent cutting properties with reduced heat generation and extremely long service life.

2009 - Lion burs, our premium diamond burs, particularly for extremely hard materials such as glass ceramics and zirconia and any other preparation.

2014 - Tungsten carbide cutters for CAD/CAM manufactured with a new generation of high precision machines.

2017 - new patents on separating strips

... Reliability

We can provide quality and innovation, since we are a family business in the 4th generation and can rely on more than 95 years of experience in the dental industry. Our success proves us right - some of our business relations have been existing for more than 80 years, and we are selling in more than 80 countries!

If you are looking for new customer-specific solutions, please contact us!

We are confident that we will convince you.

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How to use this catalogue:

Correlation between the ISO system and the HORICO® order numbers:

The ISO designation of rotary instruments consists of several groups of digits, each describing one feature.
Example:

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HORICO® order numbers for diamond instruments and carbide cutters are based on the ISO system. For diamond instruments, in the interest of simplicity letters are used for shank and grits. For carbide cutters, only the shank type is omitted. In the case of carbide burs, the old system has prevailed, which is why we use this in our catalogue, too.

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Searching systems in this catalogue:

To facilitate a quick and easy search, we offer you four search systems and a reference list:

1) Based on the item number:

Diamond instruments: The diamond instruments are sorted by the ISO shape number. If you know the HORICO order number, look for the shape number (example: FG **109**X012) and locate this ciphers under the large blue numbers above the pictures of the instruments.

Carbide cutters: The cutters are grouped according to toothing and then by shape. So you search by toothing first (e.g.: 194 **190**040) and then by shape (e.g.: **194**190 040).

2) Based on the shape:

Looking for a certain shape? For the diamond tools, there is a summary of shapes overview on page 2.

3) By type of shank:

The tables in the practice part indicate the FG or FGL-shank.

To see if there is also the head shape in FGS (shorter FG shanks), see the side reference below the table.

A collection of FGS and right-angle burs (W) can be found on the pages 39-41.

4) Based on the fields of application:

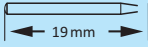
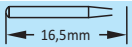
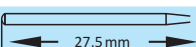
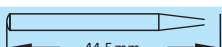
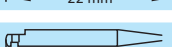
Looking for the best tool for a specific application?

For diamond burs, you find a summary overview of application recommendations from the literature on page 10. This is only meant for orientation, does not claim to be complete and does not exclude any fields of application not listed there. For carbide cutters, there is a toothing overview with application recommendations on page 76.

A reference list can be found on the inside cover. It offers the possibility of simply transposing rotary instruments from the old German designation system to the corresponding ISO designations.

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TOTAL LENGTH VARIES WITH HEAD LENGTH

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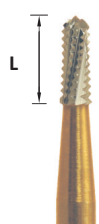
- HEADS MADE FROM FINEST TUNGSTEN CARBIDE
- HIGH PRECISION TOOTHING
- HIGHER LIFETIME THROUGH EXACT CENTRICITY
- REDUCED VIBRATIONS THROUGH STEEL SHANK

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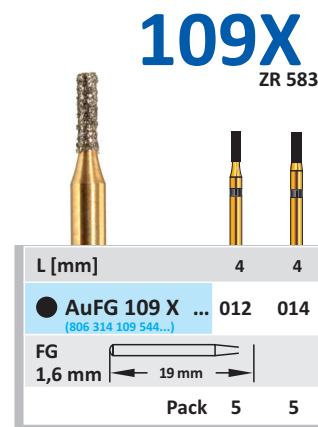
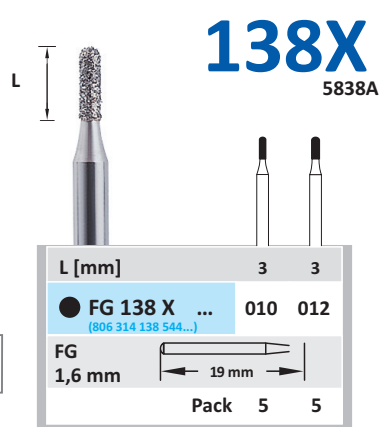
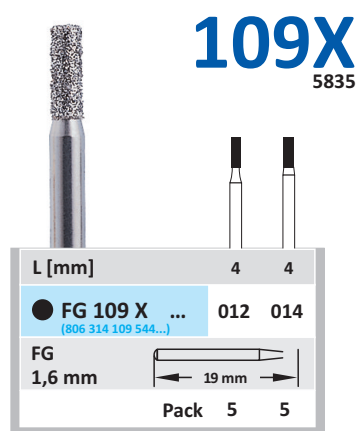
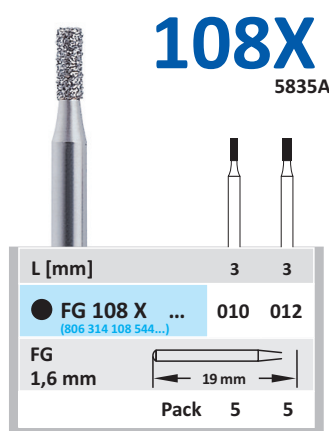
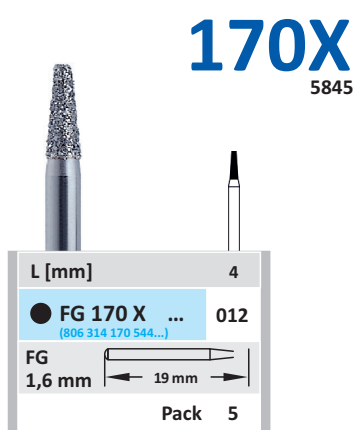
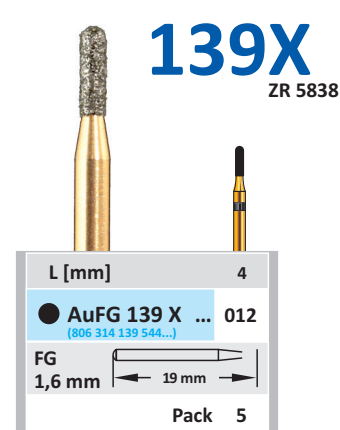
FOR CUTTING OLD FILLINGS &
CROWNS EASILY,
SHORTER CROWN CUTTER FOR MORE
SAFETY

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LION BURS

ESPECIALLY FOR ZIRCONIA AND GLASS
CERAMICS

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FOR BURS WITH DIAMETER 005 - 016

ENDODONTICS

PRIMARY MOLAR PULP THERAPY
DEFINITIVE TEETH PULP THERAPY

FIXED PROSTHODONTICS

METAL-CERAMIC CROWNS
FULL-CERAMIC CROWNS
INLAY AND ONLAY PREPARATIONS
PORCELAIN LAMINATE VENEERS

OPERATIVE DENTISTRY

COMPOSITE RESTORATIONS
CARIES DIAGNOSIS
AMALGAM RESTORATIONS
PEDIATRIC DENTISTRY
COMPOSITE RESTORATIONS FOR
PRIMARY TEETH

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INSTRUMENTS FOR
MINIMALLY
INVASIVE
PREPARATIONS

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Endodontie
traitement de la
pulpe
Odontologie
diagnostico de la
carie

ENDODONTICS
DEFINITIVE
TEETH PULP
THERAPY
**OPERATIVE
DENTISTRY**
CARIES
DIAGNOSIS

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PORCELAIN LAMINATE VENEERS

OPERATIVE DENTISTRY

COMPOSITE RESTORATIONS FOR PRIMARY

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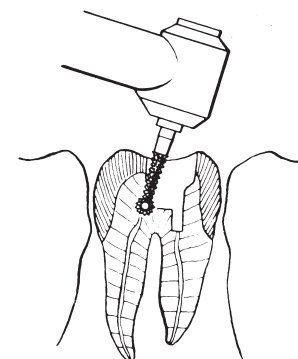
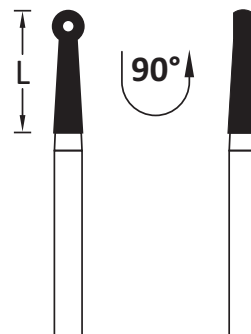
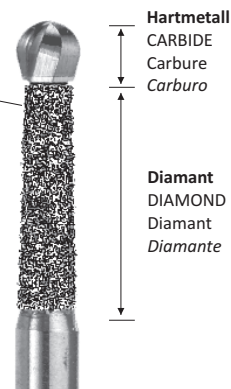
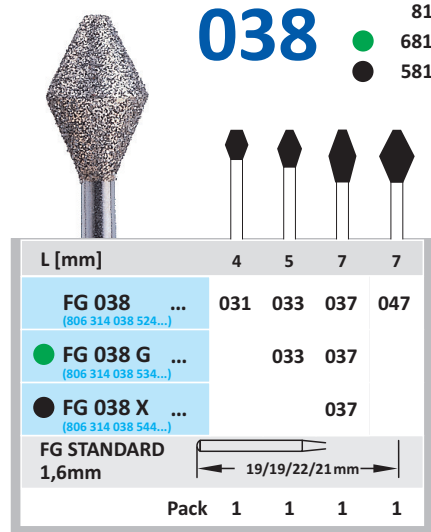
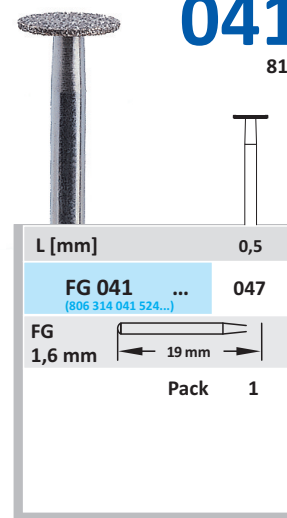
FIXED PROSTHODONTICS
PARTIAL VENEER CROWNS
METAL-CERAMIC CROWNS
ALL CERAMIC CROWNS

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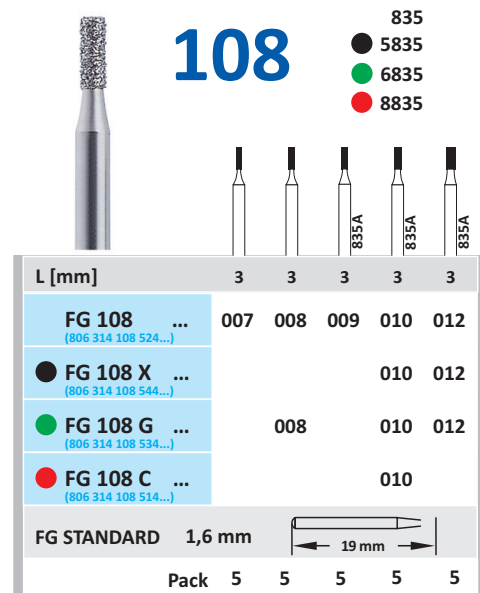
FOR CUTTING OLD FILLINGS & CROWNS
EASILY, SHORTER CROWN CUTTER FOR
MORE SAFETY

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FIXED PROSTHODONTICS
ALL CERAMIC AND METAL-
CERAMIC CROWNS

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Zahnrestauration
Komposit-Restaurationen
Amalgam-Restaurationen

OPERATIVE DENTISTRY
COMPOSITE RESTORATIONS
AMALGAM RESTORATIONS

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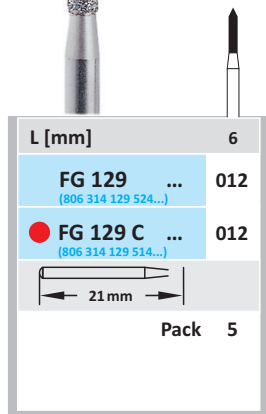
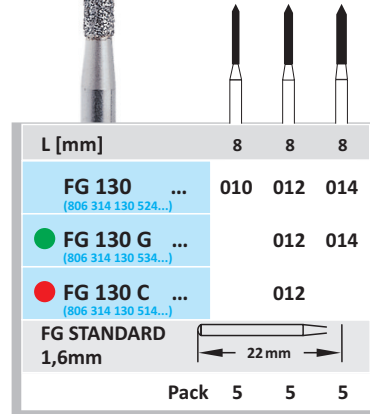
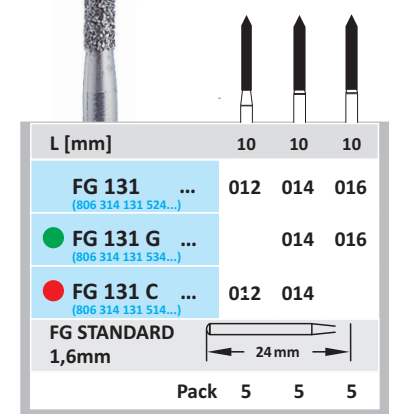
restaurations en amalgame
restaurations en composite

FIXED PROSTHODONTICS

ALL CERAMIC AND METAL-CERAMIC CROWNS
PARTIAL VENEER CROWN AND FULL VENEER CROWNS
FULL VENEER GOLD CROWN
INLAY AND ONLAY PREPARATION
METAL-CERAMIC CROWNS

OPERATIVE DENTISTRY

AMALGAM RESTORATIONS
COMPOSITE RESTORATIONS

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restauration en composite
restauration en amalgame

FIXED PROSTHODONTICS

ALL CERAMIC AND METAL-
CERAMIC CROWNS
PORCELAIN LAMINATE VENEERS
METAL RESTORATION
MOD ONLAYS
PARTIAL VENEER CROWN
OPERATIVE DENTISTRY
COMPOSITE RESTORATIONS
AMALGAM RESTORATIONS

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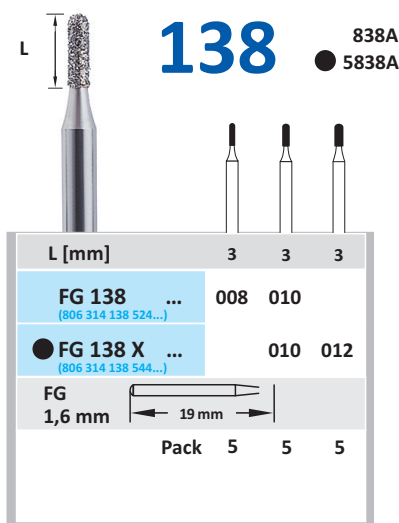
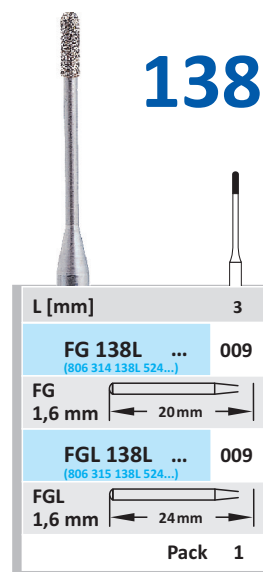
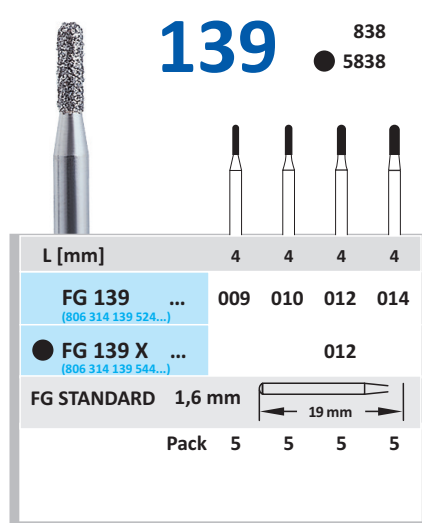
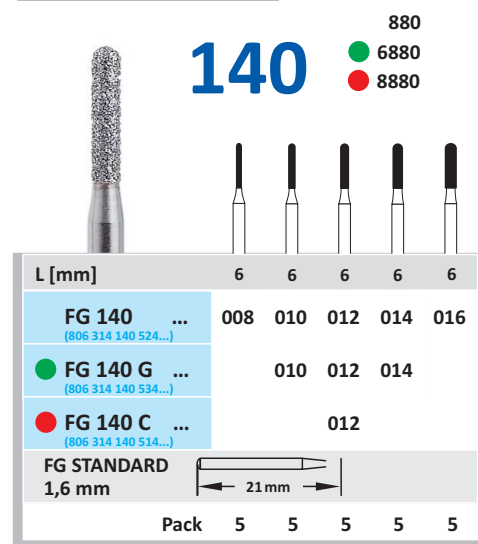
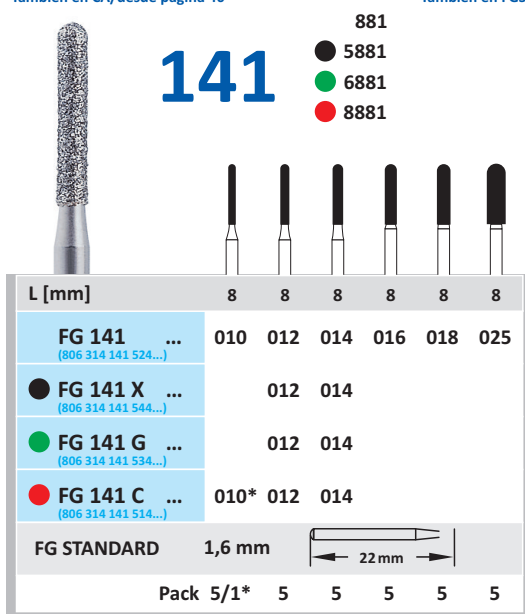
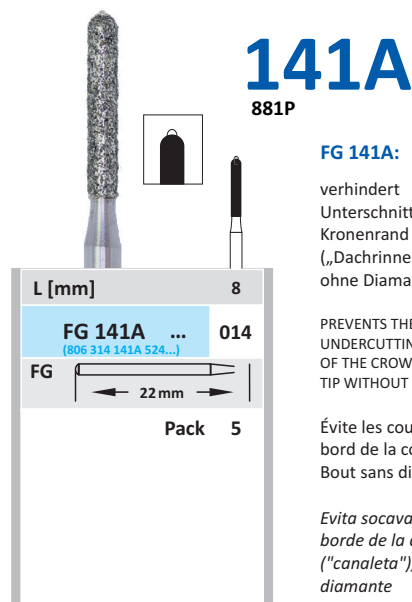
FIXED PROSTHODONTICS
ALL CERAMIC AND
METAL-CERAMIC
CROWNS

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PREVENTS THE BUR FROM
UNDERCUTTING THE EDGE
OF THE CROWN SHOULDER,
TIP WITHOUT DIAMOND

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FIXED PROSTHODONTICS
ALL CERAMIC AND METAL-
CERAMIC-CROWNS;
ADJUSTMENT OF THE
SHOULDER

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FIXED PROSTHODONTICS
ALL CERAMIC AND
METAL-CERAMIC CROWNS
PORCELAIN LAMINATE
VENEERS
METAL RESTORATION
MOD ONLAYS
PARTIAL VENEER CROWN
OPERATIVE DENTISTRY
AMALGAM RESTORATIONS
COMPOSITE RESTORATIONS

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Odontologie
restauration en composite
restauration en amalgame

FIXED PROSTHODONTICS
ALL CERAMIC AND
METAL-CERAMIC CROWNS
METAL RESTORATIONS
PARTIAL VENEER CROWNS
PORCELAIN LAMINATE VENEERS
OPERATIVE DENTISTRY
COMPOSITE RESTORATIONS
AMALGAM RESTORATIONS

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Amalgam-Restaurationen
Komposit-Restaurationen

FIXED PROSTHODONTICS

ALL CERAMIC AND
METAL-CERAMIC CROWNS
PORCELAIN LAMINATE VENEERS
METAL RESTORATIONS
MOD ONLAYS
PARTIAL VENEER CROWNS
OPERATIVE DENTISTRY
AMALGAM RESTORATIONS
COMPOSITE RESTORATIONS

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FIXED PROSTHODONTICS

PARTIAL DENTURES

FULL-CERAMIC AND METAL-CERAMIC CROWNS

PINLEDGE

INLAY AND ONLAY

OPERATIVE DENTISTRY

CARIES DIAGNOSIS

COMPOSITE RESTORATIONS

COMPOSITE RESTORATIONS FOR PRIMARY TEETH

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FIXED PROSTHODONTICS
FULL VENEER CROWNS
ALL CERAMIC CROWNS
FULL-CERAMIC
INLAY AND ONLAY
METAL-CERAMIC CROWNS
OPERATIVE DENTISTRY
COMPOSITE RESTORATIONS

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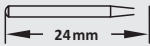
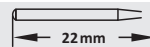
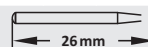
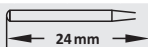
ORAL AND MAXILLOFACIAL SURGERY
TOOTH EXTRACTION

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FIXED PROSTHODONTICS
ALL CERAMIC AND METAL-CERAMIC
OPERATIVE DENTISTRY
COMPOSITE RESTORATIONS

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OPERATIVE DENTISTRY
AMALGAM RESTORATIONS
COMPOSITE RESTORATIONS
ALSO FOR PRIMARY TEETH

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FOR CUTTING OLD FILLINGS &
CROWNS EASILY, SHORTER CROWN
CUTTER FOR MORE SAFETY

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FIXED PROSTHODONTICS
CUSTOM CAST DOWEL CORES
FULL-CERAMIC CROWNS
METAL CERAMIC CROWNS
MOD ONLAYS
PORCELAIN VENEERS
PROXIMO-OCCLUSAL INLAYS
OPERATIVE DENTISTRY
CARIES DIAGNOSIS
COMPOSITE RESTORATIONS

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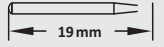
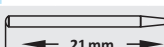
PREVENTS THE BUR FROM UNDERCUTTING THE EDGE OF THE CROWN SHOULDER, TIP WITHOUT DIAMOND

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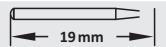
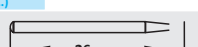
FIXED PROSTHODONTICS
FULL-CERAMIC CROWNS
METAL CERAMIC CROWNS
PORCELAIN VENEERS
PROXIMO-OCCLUSAL INLAIS
OPERATIVE DENTISTRY
CARIES DIAGNOSIS
COMPOSITE RESTORATIONS

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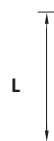
FIXED PROSTHODONTICS
ALL CERAMIC CROWNS
METAL CERAMIC CROWNS
PARTIAL VENEER CROWN

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FIXED PROSTHODONTICS
METAL CERAMIC CROWNS
OPERATIVE DENTISTRY
COMPOSITE RESTORATIONS

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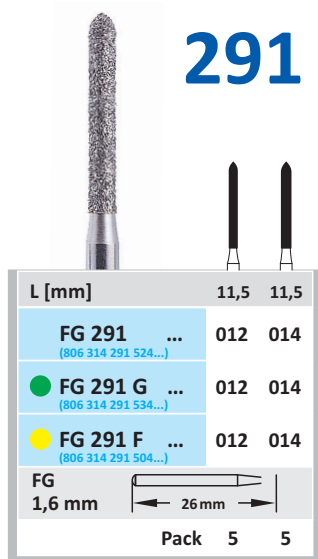
FIXED PROSTHODONTICS
ALL CERAMIC AND METAL-CERAMIC
ALL METAL RESTORATIONS
CERAMIC LAMINATE VENEERS
METAL-CERAMIC RESTORATIONS
OPERATIVE DENTISTRY
AMALGAM RESTORATIONS

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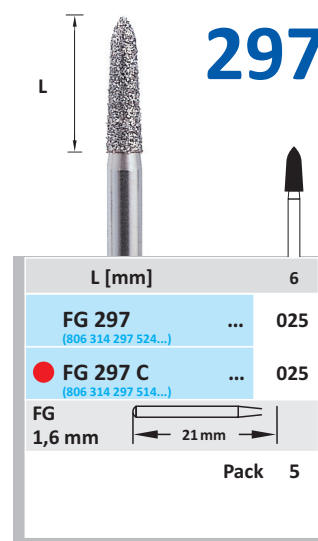
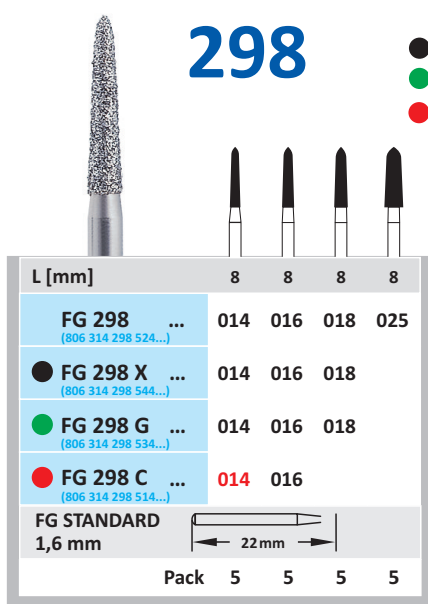
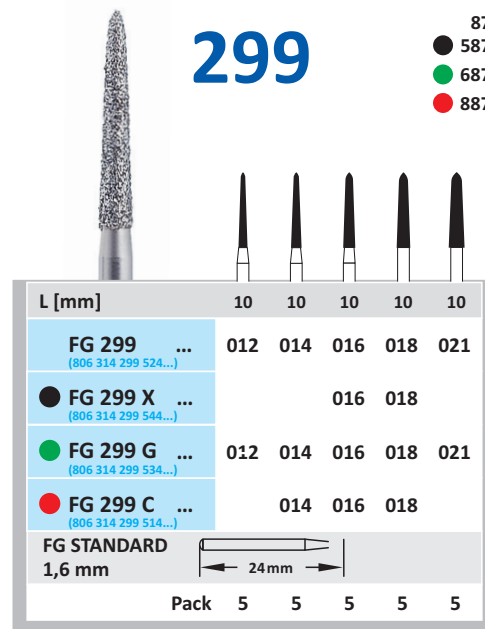
FIXED PROSTHODONTICS
METAL-CERAMIC RESTORATIONS
ALL CERAMIC AND METAL-CERAMIC
PARTIAL VENEER CROWNS

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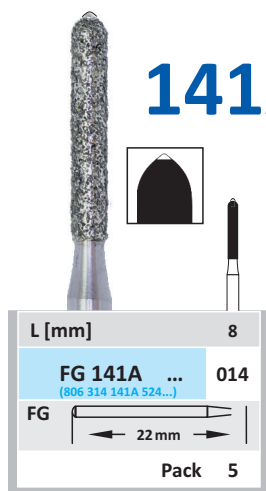
PREVENTS THE BUR FROM
UNDERCUTTING THE EDGE
OF THE CROWN SHOULDER

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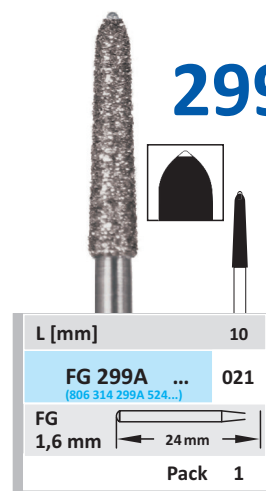
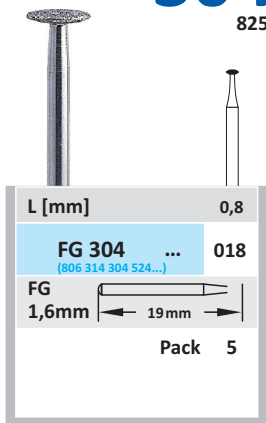
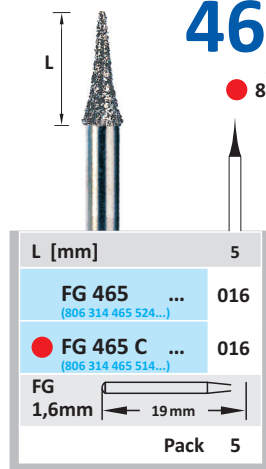
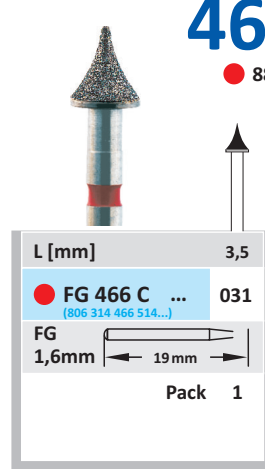
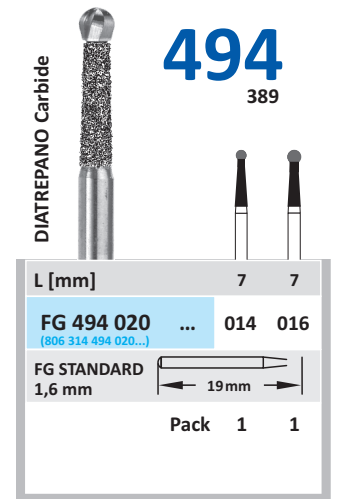
FG 465 + FG 466:

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FIXED PROSTHODONTICS
FOR FISSURES

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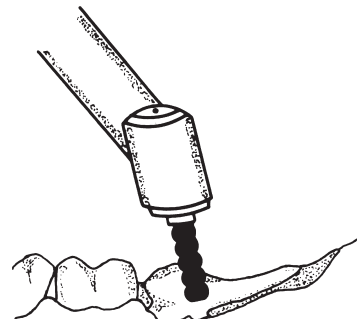
ORAL AND MAXILLOFACIAL
SURGERY
REMOVAL OF MISALIGNED
WISDOM TEETH

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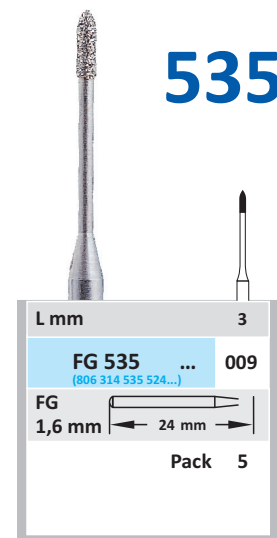
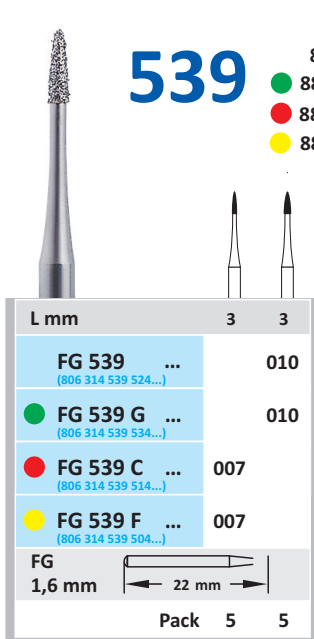
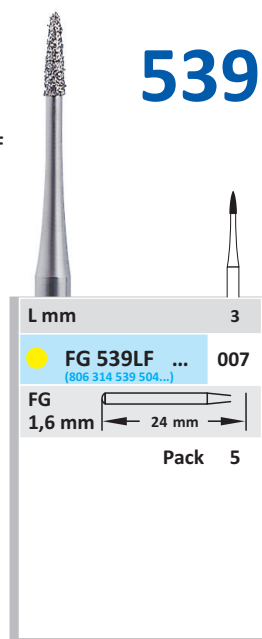
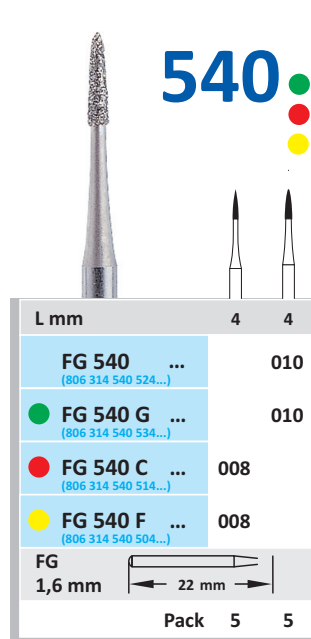
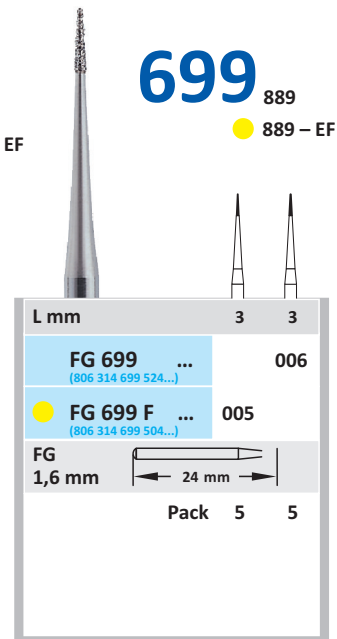
FIXED PROSTHODONTICS

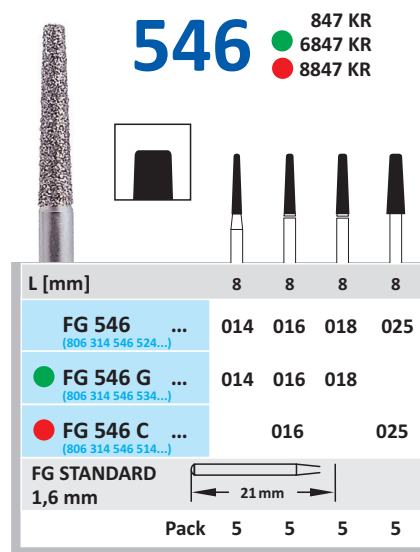
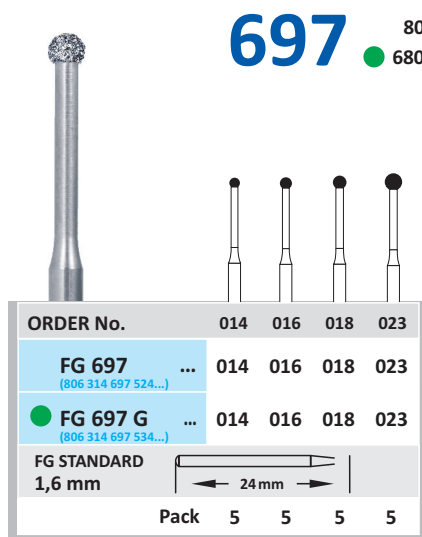
FULL VENEER CROWNS
ALL CERAMIC
FULL-CERAMIC CROWNS
INLAY AND ONLAY
METAL-CERAMIC CROWNS
OPERATIVE DENTISTRY
COMPOSITE
RESTORATIONS

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ENDODONTICS

DEFINITIVE TEETH PULP THERAPY
OPERATIVE DENTISTRY
CARIES DIAGNOSIS

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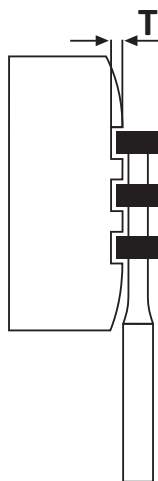
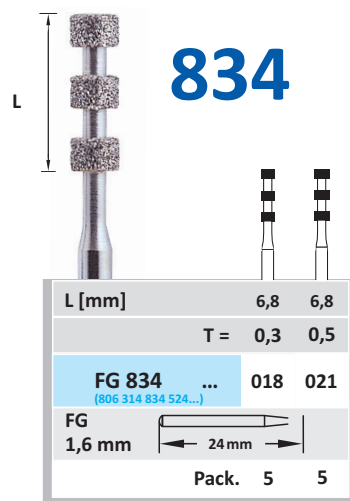
(PRAXIS) DIAMANTSCHLEIFER - ORALCHIRURGIE

(SURGERY) DIAMOND BURS - OPERATIVE DENTISTRY

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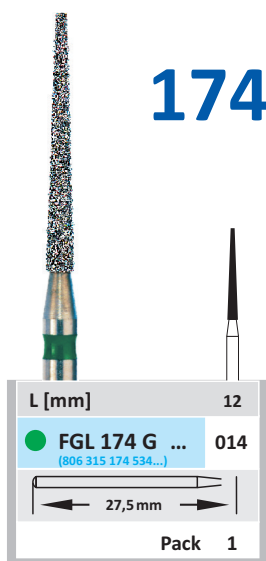
ORAL AND MAXILLOFACIAL
SURGERY
TOOTH EXTRACTION

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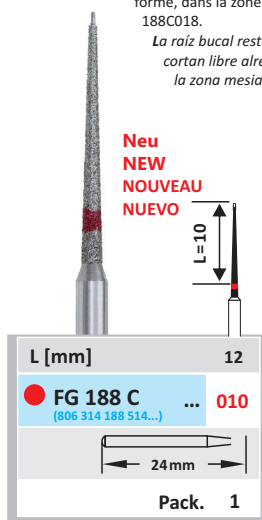
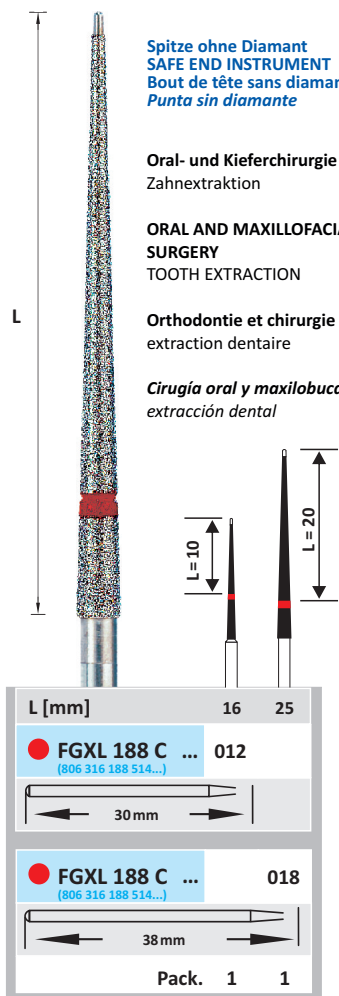
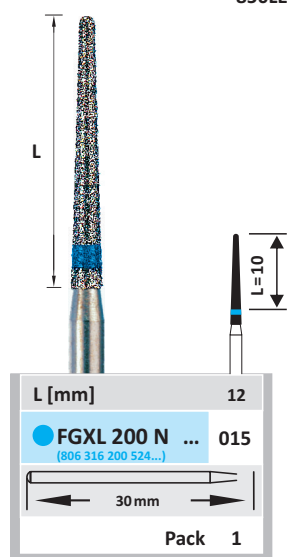
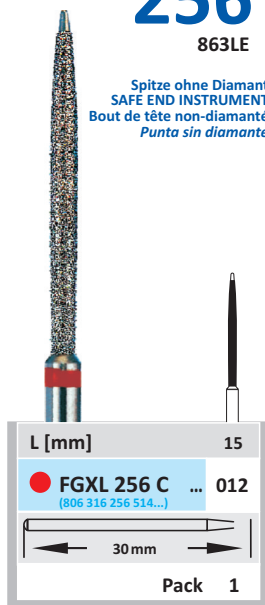
ORAL AND
MAXILLOFACIAL
SURGERY
TOOTH EXTRACTION

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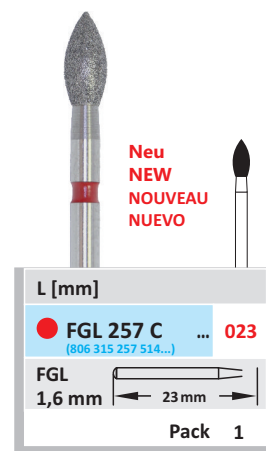
FIXED PROSTHODONTICS
METAL-CERAMIC
RESTORATIONS
ORAL AND
MAXILLOFACIAL
SURGERY
TOOTH EXTRACTION

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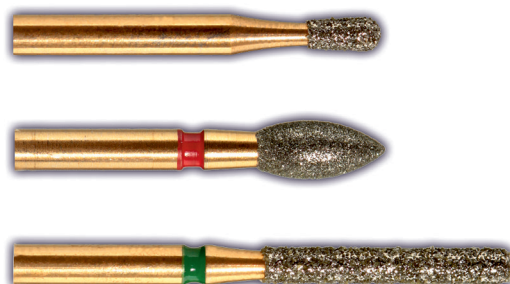
DESIGNED FOR ZIRCONIA AND LITHIUM-(DI)-
SILICATE, BUT ALSO OUTSTANDING ON ALL
OTHER DENTAL PREPARATIONS.

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YOUR BIG ADVANTAGES COMPARED TO STANDARD-DIAMONDS:

- THEY WORK
- FAST AND DURABLE ON HARD CERAMICS LIKE ZIRCONIA AND LITHIUM-DI-SILICATE
- CUT EVEN FASTER – REDUCED CHAIRTIME THROUGH CUTTING EFFICIENCY, WHICH IS INCREASED BY AT LEAST 20% (MM/MIN) ON ALL OTHER MATERIALS
- COOLER COMPARED TO STANDARD BURS IN TREATMENTS WITH THE SAME DURATION
- EXCELLENT PRICE/PERFORMANCE RATIO

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OUR TURBO DIAMONDS (TFG) OF THE LION FAMILY WORK:

- **FASTER**, BECAUSE THE RIM OF EACH SPIRAL IS WORKING LIKE AN ADDITIONAL BLADE.
- **EASY AND PERSISTENT**, BECAUSE THE ABRASIVE WEAR IS TRANSPORTED FASTER IN GOLDEN GOUGES AND WITHOUT CLOGGING THE DIAMOND GRIT.
- **WITH AN IMPROVED SAFETY**, BECAUSE THE SPIRALS ARE POSITIONED IN A WAY THAT THE SHANK OF THE BUR GETS PRESSED INTO THE TURBINE.

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FOR DIFFICULT TO REACH AREAS AND PEDIATRIC DENTISTRY

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DIAMOND BURS WITH RIGHT ANGLE AND FG SHORT SHANK

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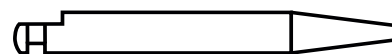
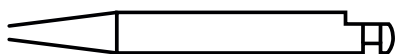
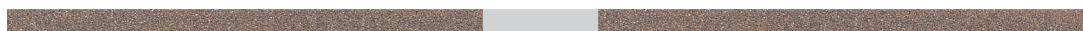
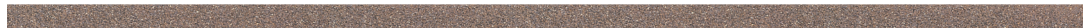
THESE STRIPS HAVE A MUCH LONGER LIFETIME
 BECAUSE OF THEIR DIAMOND COATING COMPARED
 TO STEELCARBO®- STRIPS WITH A RUBY LAYER.
 THEY ARE AUTOCLAVABLE AND STERILIZABLE.

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ALL OUR STRIPS ARE BASED ON STAINLESS STEEL
 AND THUS AUTOCLAVABLE AND STERILIZABLE.

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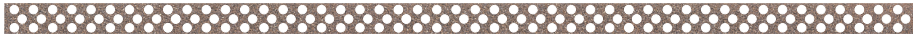
STRIPS HAVE A MUCH LONGER LIFETIME DUE TO THEIR DIAMOND COATING COMPARED TO STEELCARBO® STRIPS WITH A RUBY LAYER. TRANSPORT OF DEBREE IS BETTER BECAUSE OF THE PERFORATION. THEY PREVENT TOOTH SUBSTANCE MUCH BETTER FROM DAMAGE BY HEAT THROUGH THE COOLING EFFECT. THEY ARE AUTOCLAVABLE AND STERILIZABLE.

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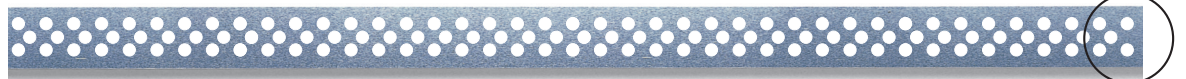
IN ADDITION TO THE ABOVE WRITTEN ADVANTAGES THE GAP ENABLES AN EASIER "START" OF THE PREPARATION. ALL STRIPS ARE MADE OF STAINLESS STEEL AND CAN BE STERILIZED AND AUTOCLAVED.

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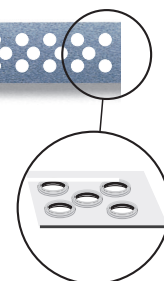
ABRASIVE STAINLESS STEEL STRIPS
Tiras de profilaxis

VERY THIN STRIP MADE FROM STAINLESS STEEL FOR PROPHYLAXIS

- 1) GENTLE REMOVAL OF PLAQUE, STAINS AND RESIDUES OF FILLINGS
 - 2) DO NOT CONTAIN DIAMOND GRIT. THE STEEL EDGES AROUND THE PERFORATION CUT WITHOUT HARMING THE ENAMEL
 - 3) THE PERFORATION PICKS UP ABRASION AND EXCESS MATERIAL
 - 4) IS MUCH MORE TEAR-RESISTENT AND DURABLE THAN THE NORMALLY USED STRIPS MADE OF PLASTICS.
- AUTOCLAVABLE AND STERILIZABLE.

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GENUINE CORUNDUM GALVANICALLY COATED ON STAINLESS STEEL,
CAN BE CUT TO DESIRED LENGTH BEFORE USE.
AUTOCLAVABLE AND STERILIZABLE.

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GENUINE CORUNDUM GALVANICALLY COATED ON STAINLESS STEEL,
WITH INSERTION GAP.
AUTOCLAVABLE AND STERILIZABLE.

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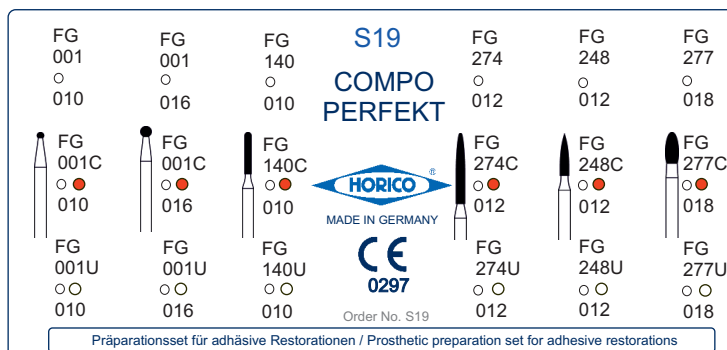
Operative preparation set for adhesive restorations

Indication:

- Cavity preparation procedures and shaping of composite filling material
- Anterior composite filling restorations
 - Anterior incisal edge and cusp build-ups
 - Adhesive bridges (Maryland type)
 - Posterior ceramic inlays/onlays
 - Aesthetic veneer preparations for either ceramic or composite material

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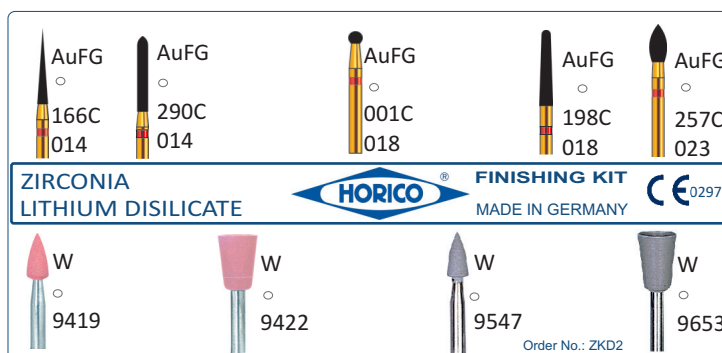
Finishing-Set for crown & bridges made from zirconia and glass ceramics

Indication:

- to refine and adjust crown margins and interproximal areas
- to adapt occlusal contact areas
- to refine and polish restoration surfaces

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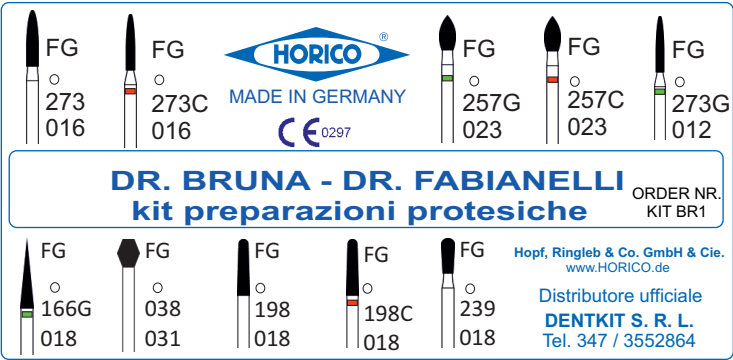
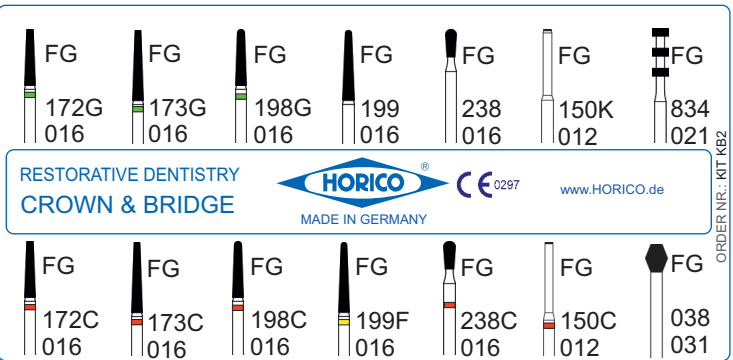
KIT FOR PROSTHETIC PREPARATION
PARTICULARLY SUITABLE FOR THE
VERTICAL TECHNIQUE OF
CURETTAGE.

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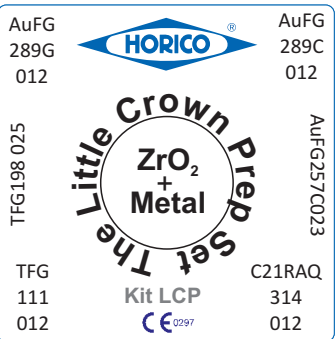
KIT FOR DENTAL PROSTHETICS,
IDEAL FOR ALL CROWN AND
BRIDGE PREPARATIONS.

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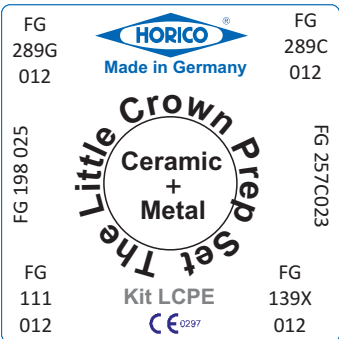
KITS FOR CROWN PREPARATION.
LCP: WITH TFG 198 025 FOR FAST REDUCTION AND
LION BURS, WHICH CAN BE ALSO USED FOR ZIRCONIA.
LCPE: FG 139X012 FOR REMOVAL OF OLD CROWNS.

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- ONLY ONE KIT

- ONLY ONE KIT FOR ALL DENTURE ADJUSTMENTS
- THE SPECIAL BLACK COATING LEADS TO A LONG LIFETIME OF TC-CUTTERS AND REDUCES CLOGGING

RATING:

EVALUATION BY THE
INDEPENDENT US REVIEW
„CLINICIANS REPORT®“: 100%
PARTICIPATING DENTIST
WOULD RECOMMEND THE KIT
DA TO COLLEAGUES. REPORTS
CAN BE ORDERED AT HORICO

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Use on all other dentures

S274
190
060
○

**Denture Adjustment Kit
for all materials**

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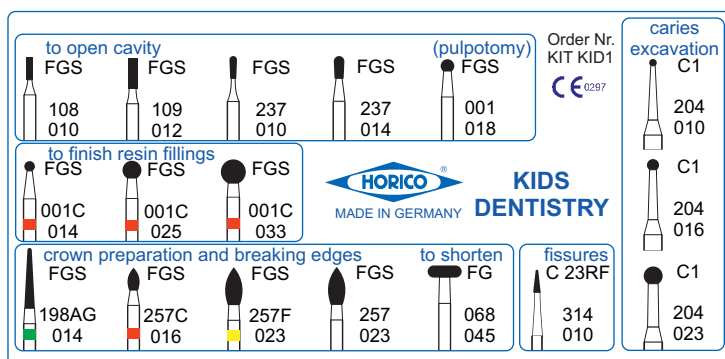
2) CROWN PREPARATION (OPEN UP THE APPROXIMAL SPACE: FGS198AG014; SHORTENING OF OCCLUSAL AREAS: FG068 045; CHAMFER THE EDGES: FGS257C016, FGS257F023, FGS257 023).

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FOR FURTHER
INFORMATION,
PLEASE REFER
TO OUR
CATALOGUES
FOR POLISHERS
AND TC-BURS.

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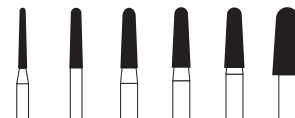
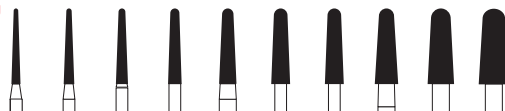
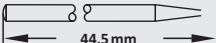
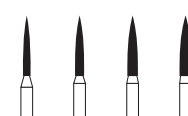
FG FOR LABORATORY TURBINE - WITH WATER COOLING

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DESIGNED FOR ZIRCONIA AND
GLAS CERAMICS

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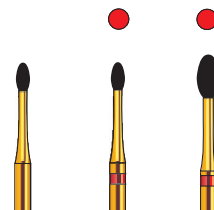
2° FG-DIAMOND INSTRUMENTS
FOR THE PRECISION MILLING
TECHNIQUE

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MATERIALS: CERAMICS, VENEER
ACRYLICS, POLYURETHAN, RESIN,
SOFT ACRYLICS

Matériaux:
polyuréthane, céramique
polyuréthane, résines acryliques,
résines molles
Material:
cerámica, blindajes de plástico,
poliuretano, plásticos, rebases
blandos

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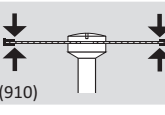
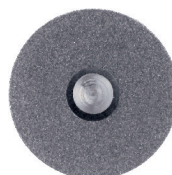
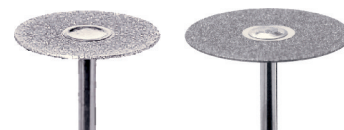
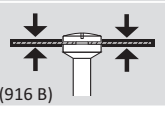
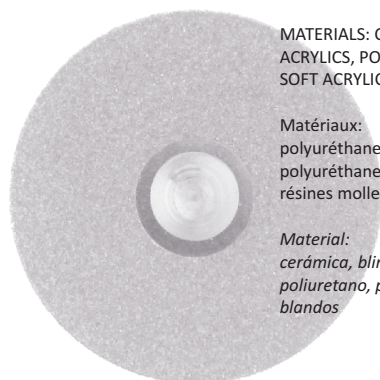
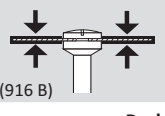
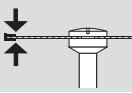
PLASTER SECTIONING DISCS FOR MODEL PREPARATION, TEETH AVOID CLOGGING AND ALLOW COOL WORKING

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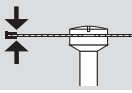
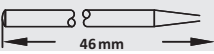
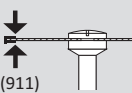
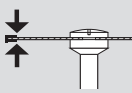
MATERIALS: CERAMICS, PLASTER, POLYURETHAN, RESIN, SOFT ACRYLICS, VENEER ACRYLICS

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PLASTER SECTIONING DISCS FOR MODEL PREPARATION

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PLASTER SECTIONING DISCS
FOR MODEL PREPARATION

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MATERIALS: CERAMICS,
PLASTER, POLYURETHAN,
RESIN, SOFT ACRYLICS,
VENEER ACRYLICS

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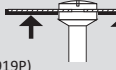
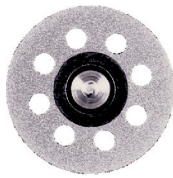
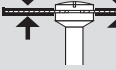
FOR TREATING
ZIRCONIA WITH
WATER COOLING.

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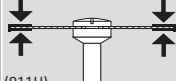
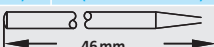
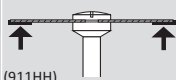
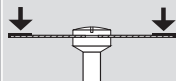
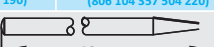
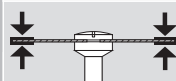
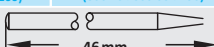
H 355G:

Material: Keramik, Gips,
Polyurethan, Kunststoffe,
Verblendkunststoffe

MATERIALS: CERAMICS,
PLASTER, POLYURETHAN,
RESIN, VENEER ACRYLICS

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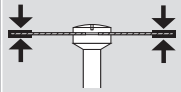
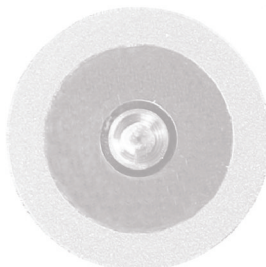
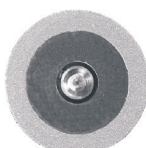
FOR CUTTING OFF SPRUES FROM PRESSED CERAMICS

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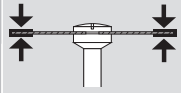
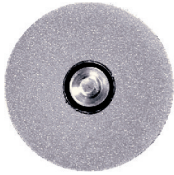
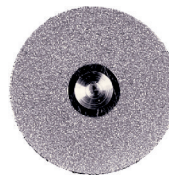
SERRATED PLASTER
SECTIONING DISC FOR MODEL
PREPARATION, AVOIDS
CLOGGING.

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MATERIALS: CERAMICS,
PLASTER, POLYURETHANE,
RESIN, SOFT ACRYLICS, VENEER
ACRYLICS

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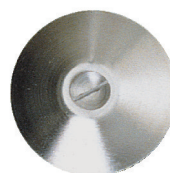
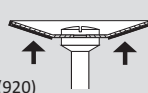
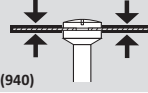
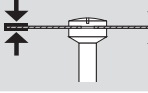
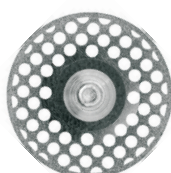
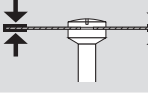
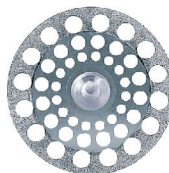
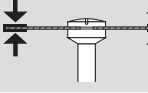
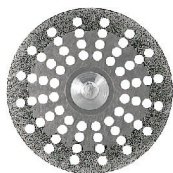
H 934 G:

Material: Keramik, Gips, Polyurethan,
Kunststoffe, Kunststoffverblendungen

MATERIALS: CERAMICS, PLASTER,
POLYURETHAN, RESIN, SOFT ACRYLICS,
VENEER ACRYLICS

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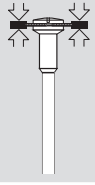
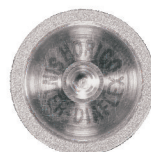
MATERIALS: CERAMICS, GOLD,
STEEL, ALL CHROME-COBALT
ALLOYS

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INSTRUCTIONS FOR USE:

TO AVOID BREAKAGE PLEASE USE
WITHOUT PRESSURE.
ALWAYS WEAR PROTECTIVE GLASSES!
DO NOT USE INSIDE MOUTH!
TO BE SHARPENED AND CLEANED WITH DRY
DIACLEAN (ORDER NO. 500).
DO NOT REMOVE FROM MANDREL!

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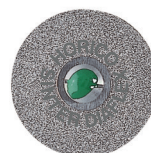
ATTENTION! MARKED YELLOW
FOR SINGLE SIDED DISC!

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RAPID CLEANING OF CLOGGED DIAMONDS WITH
DIACLEAN. MOISTEN DIACLEAN WITH WATER AND
GRIND THE CLOGGED DIAMOND INSTRUMENT INTO IT.

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INSTRUCTIONS FOR USE:

TO AVOID BREAKAGE PLEASE USE
 WITHOUT PRESSURE.
 ALWAYS WEAR PROTECTIVE GLASSES!
 DO NOT USE INSIDE MOUTH!
 TO BE SHARPENED AND CLEANED WITH DRY
 DIACLEAN (ORDER NO. 500).
 DO NOT REMOVE FROM MANDREL!

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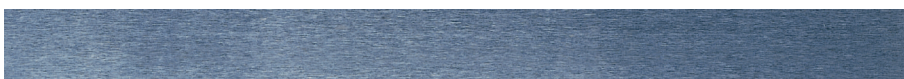
**MATERIALS: METAL ALLOYS, CERAMICS,
 VENEER ACRYLICS**

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SINTERED BURS FOR THE LABORATORY TURBINE

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- DIAMOND, SYNTHETIC BONDING
- HIGH PERFORMANCE
- MINIMAL HEAT DEVELOPMENT, THEREFORE NO WATERCOOLING REQUIRED

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FOR CERAMICS, METAL CERAMICS AND METALS

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- DIAMOND, CERAMIC BONDING
- HIGH SHAPE STABILITY
- HEAT RESISTANT
- DURABLE, FINE GRIT

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RAPID CLEANING OF CLOGGED DIAMONDS WITH DIACLEAN. MOISTEN DIACLEAN WITH WATER AND GRIND THE CLOGGED DIAMOND INSTRUMENT INTO IT.

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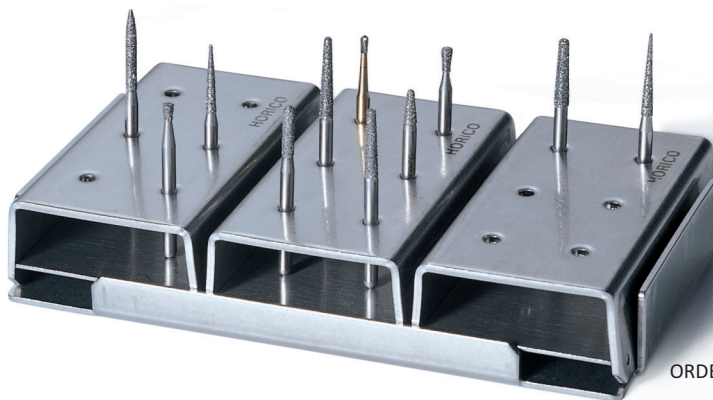
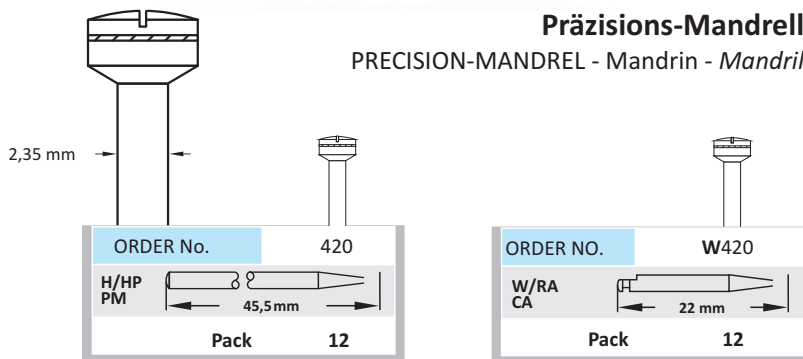
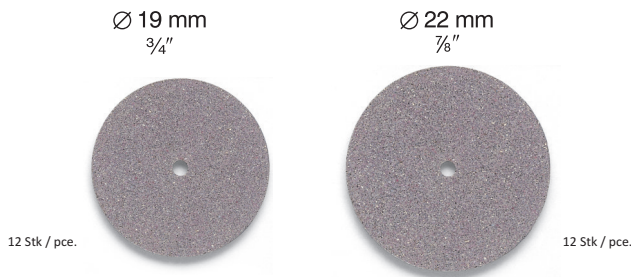
GENUINE CORUNDUM GALVANICALLY BONDED ON STAINLESS STEEL. THESE DISCS ARE AUTOCLAVABLE AND STERILIZABLE.

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To process non precious metals from master model
to polished model cast frame.

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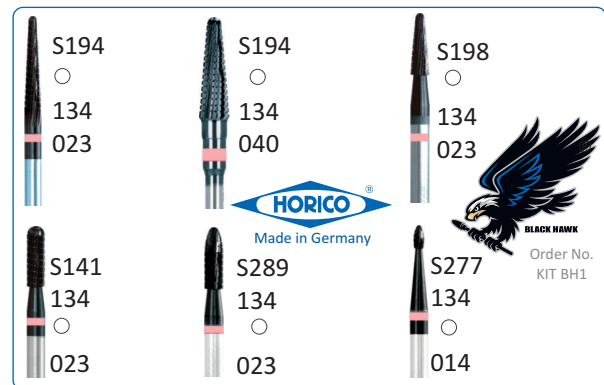
Selection of Black Hawk Cutters with toothing 134 to
process non precious metal alloys like cobalt-chrome,
nickel-chrome or titanium.

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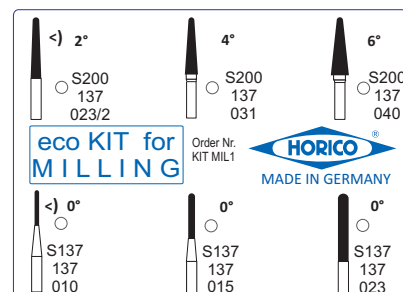
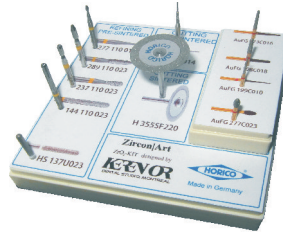
Selection of Black Hawk Cutters with toothing 137 for the
milling technique in 0°, 2°, 4° und 6°.

Applicable for non precious metal alloys like chromium-
cobalt, chromium-nickel or titanium.

Zirkonset Zircon | Art (SIKER2)

Designed by KerenOr Dental Studio, Montreal

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Scope of this information

Unless stated otherwise, the following information applies to the groups of products below offered by Hopf, Ringleb & Co. GmbH & Cie. under the brand name HORICO® when used as intended: Separating strips, Diamond burs, Carbide burs, Polishers

Intended use

All the instrument groups listed above are for multiple use during dental treatments. They are used for reduction, excavation and cutting or for surface conditioning of dental tissues such as bone, dentin and dental enamel and dental materials such as composites, metal alloys, ceramics and dental resins.

They are intended for use by trained technicians and dentists. Separating strips are used manually. All other instruments are powered by a handpiece or angle attachment.

Note: There are handpiece and angle attachment drives with collets for handpiece, right angle shank and friction grip (FG) shanks. These must comply with the relevant ISO standards.

Transport, storage and treatment of new instruments prior to first use

There are no particular conditions necessary for transport.

Please store all original packaged instruments at a clean and dry place at room temperature.

New instruments are not packaged in sterile condition and must be cleaned, disinfected and sterilized according to preparation instructions before each use.

With unpacking longer, thinner instruments and discs, take care not to bend them. It is best to open blister packaging with scissors; with rigid packages, both labels should be cut before lifting the cap.

Disposal

All our instruments are free of hazardous substances ex works. Only the requirements for protection against infection and contamination need be considered for disposal.

Instructions for use

All our instruments are developed and manufactured with care for their particular application. Improper use can endanger the user, the patient and possibly other persons as well as damage the instrument and handpiece.

Prior to use, ensure that:

- Users and assistance where mouth protection, goggles and gloves and that the environment (treatment unit, etc.) is appropriately disinfected, because infection particles can be scattered by the fast rotation and spray water.
- Use only drives which are in good technical and hygienic condition. Please follow the operating instructions from the handpiece manufacturer!
Please note explicitly that most preprocessing units generally do not lubricate the collet and bearings. Most handpiece manufacturers usually require this after 20 to 30 minutes of operation.
- Please clamp the shaft of the instrument as deeply as possible. If instruments are loose or protrude too far, they can fly off, bend or break, causing injuries, or be swallowed or aspirated.
- Select the speed such that the maximum allowed RPM is not exceeded (see the table for maximum speed). Exceeding the maximum speed increases safety risks, reduces the quality of work and generates heat. The recommended speed, which is generally about half the maximum speed, produces the best work results and reduces undesirable secondary effects to a minimum.
- Ensure sufficient air/water cooling (minimum 60 ml/min).
- Processing extraneous materials in the mouth, such as filling materials, can release nanoparticles of these. Depending on the starting material, these may be bioactive. Thus suction and possibly other protective measures such as a dental rubber dam, etc. are recommended.
- Please bring the instrument up to working speed outside of the mouth or prior to contact with the workpiece. If vibrations occur, the instrument is bent and can no longer be used!
- Please work with as little pressure as possible (about 50 g, corresponding to the pressure applied when writing) and do not twist the instrument. Higher pressure only leads to greater heat development, faster wear and inferior work results.
- If the instrument jams, stop the handpiece, carefully remove the instrument without twisting it and check for damage before using it again.
- Dull and damaged instruments must no longer be used. Please check the instruments prior to each use.
Signs of damage with diamond instruments are blank spots on the working part, bent instruments which produce vibrations when starting and changes in the original form. Carbide instruments exhibit damaged and deformed cutting edges or breaks.
- Long instruments are not suited for canals with curvature: there is a risk of breakage.

If the safety instructions are not followed, damage can occur to the tooth and surrounding tissue or the workpiece, possibly endangering the user, the patient and other persons.

Maximum speed
Diamant-Schleifinstrumente - DIAMOND INSTRUMENTS
INSTRUMENTS DIAMANTES - INSTRUMENTOS DIAMANTADOS

ISO Ø [$\frac{1}{10}$ mm]	- FG, W, RA, CA - [UPM - RPM - TPM]	- H, HP, PM - [UPM - RPM - TPM]
005-014	450 000	50 000
016-023	300 000	50 000
025-045	120 000	50 000
047-065	80 000	50 000
066-093	60 000	40 000
100-127	30 000	30 000
130-300	25 000	25 000
450		20 000

Sinter-Diamanten -
SINTERED DIAMONDS
DIAMANTS FRITTES
DIAMANTADOS SINTERIZADOS

ISO Ø [$\frac{1}{10}$ mm]	- H, HP, PM - FG - [UPM - RPM - TPM]
016-070	25 000
080-310	15 000

Hartmetall-Instrumente
TUNGSTEN CARBIDE INSTRUMENTS
INSTRUMENTS EN CARBURE DE TUNGSTENE
INSTRUMENTOS DE CARBURO DE TUNGSTENO

ISO Ø [$\frac{1}{10}$ mm]	- FG - W, RA, CA - [UPM/RPM/TPM]	- H, HP, PM - [UPM/RPM/TPM]
003-027	200 000	50 000
031	120 000	50 000
040	70 000	50 000
045	65 000	50 000
050	60 000	40 000
060	50 000	35 000
070	30 000	30 000

The optimal working speed recommended by HORICO is usually about 50% of the maximum permitted speed. However, the speed should be adjusted according to the material (refer to the information for the item in the catalog).

Note regarding the obligation to report:

According to the relevant statutory regulations, serious incidents which occur in conjunction with a product from us must be reported to us (at the address below) and to the responsible national authorities of the country in which the user works.

In terms of the law, incidents are subject to mandatory reporting if the use of the instrument directly or indirectly caused or could have caused the death or serious damage to health of the patient, user or another person (source: the German Safety Plan for Medical Devices [MPSV], Section 2).

Reprocessing instructions

General information

The following instructions contain the minimum requirements for reprocessing of which have been validated for safety. In addition, there are processing steps marked as recommendations (in italics), which can improve the result of the controls, particularly after cleaning and disinfection. Furthermore, there is also a validated manual method provided for cleaning and disinfection. It is considered as an alternative for regions with poor infrastructure or for cases where technical equipment fails. Fundamentally, all equipment used must meet the relevant standards and regulations and be in proper technical and hygienic condition.

Basic legal and regulatory information:

These preparation instructions were produced in accordance with the EU medical device regulation (MDR/2017), DIN EN ISO 17664: 2004 and the recommendations of the commission for hospital hygiene and infection prevention at the Robert Koch Institute (RKI): "Infection prevention in dentistry – requirements for hygiene" (2006) and "Hygiene Requirements for the Reprocessing of Medical Devices" (2012)

Scope of this information

Unless stated otherwise, the following information applies to the groups of products below offered by us when used as intended: Separating strips, Diamond burs, carbide burs, Polishers.

Risk assessment and classification of the product groups according to the guideline cited above

Due to their area of application and abrasive surface and/or blind holes, separating strips and polishers are to be classified as semi-critical B for polishing, with particular requirements for preparation. Diamond burs and carbide burs are to be classified as semi-critical B (nonsurgical) or critical B (surgical) depending on the application and due to services with particular requirements for preparation. All the aforementioned product groups are technically suitable for preparation in most types of ultrasonic baths, cleaning and disinfection equipment and autoclaves. Cleaning and disinfection equipment can be insufficient for sole use in cleaning under some circumstances due to the aforementioned special considerations. (Please refer to critical process steps as well.) Our instruments are not suited for Chemiclav or hot air sterilization.

Critical process steps

The post-cleaning check is to be considered a critical process step (critical control point). The rough surfaces from diamond coating or serration of instruments and blind holes in polishers result in a risk that tissue remains and other impurities are not removed completely. This can reduce the effectiveness of disinfection and sterilization. Therefore, drying of impurities is to be avoided, and the post-cleaning check should be performed with particular care using aids such as a magnifier or microscope.

Expertise and protection of the person performing preparation

The aforementioned product groups must only be prepared by persons with the necessary expertise. The requirements for this are defined in the relevant regulations and statutes. Gloves, mouth protection, safety glasses and protective clothing must be worn to avoid contact with potentially infectious material.

Suitable media for cleaning and disinfection

Only certified media produced according to the relevant standards may be used. The manufacturers of these cleaning and disinfection materials check whether a product is fundamentally usable for rotary instruments. Thus all media for which the manufacturer indicates a corresponding approval can be used. They should contain corrosion protection agents. Cleaning and disinfection media which contain strong acids and bases or hydrogen peroxide and sodium hypochlorite, are not suited for use.

Storage and preparation of new instruments

Please store all original packaged instruments in a clean and dry place at room temperature. New instruments are not packaged in a sterile condition and must be disinfected and sterilized according to the following instructions before use.

Preparation of used instruments

All product groups must be prepared prior to each use according to the following steps.

1. Cleaning and disinfection

Cleaning should begin as soon as possible after use, no later than 6 hours, to avoid drying contaminants which are difficult to remove. Before storing used instruments, procedures for occupational safety (involving the risk of infection and contamination) must be followed.

Recommendation: All types of burs should be stored in a "bur's bath" as soon as possible after use until reprocessing. The bur bath should contain a disinfectant for better protection of the person responsible for reprocessing.

A) Automated cleaning and disinfection:

Recommendation: Pre-cleaning in an ultrasonic bath can significantly improve visual inspection results. Please note the information from the device manufacturer.

Clean and disinfect the cleaning and disinfection equipment according to the instructions of the manufacturer. It should be ensured that the instruments are held in place so that they do not rub against anything and in particular so that the heads are not in a dead zone for rinsing. Instruments with blind holes (polishing equipment) should be stored flat.

B) Manual cleaning and disinfection

Place instruments in cold water for at least five minutes, rinse under flowing water, remove any tissue remains with the aid of a brush or DIACLEAN cleaning stone (the latter for diamond grinding instruments only), then rinse again. Repeat brushing and rinsing as needed after inspection with a magnifier.

After that, disinfect strip, diamond and carbide instruments in a chemical disinfection bath. Please use only disinfectants which contain corrosion protection agents and comply with the concentrations and treatment times prescribed by the manufacturer.

All instruments must be rinsed after disinfection with water which is microbiologically and chemically safe, and then dry thoroughly with medical-grade compressed air.

2. Follow-up inspection

Prior to sterilization, all instruments must be inspected for cleanliness and functionality, preferably with a magnifier or microscope. For each kind of impurity, particularly on the diamond coating or serration or in blind holes (for polishers), the cleaning and disinfection process must be repeated or the instrument must be discarded. Worn, bent or damaged instruments must be excluded.

3. Sterilization

Recommendation for semi-critical (nonsurgical) instruments used: we recommend unpackaged sterilization in autoclaves as a minimum.

Critical (surgical) instruments used must be packaged according to regulations and marked for sterilization in the autoclave. Generally, a partial pre-vacuum and dwell time of five minutes at 134°C and about 2 bar pressure are used. The information from the equipment manufacturer is authoritative for this.

Note: Galvanic corrosion can occur in autoclaves with steel holders. This can lead to surface rust on the stainless steel shafts of instruments which, however, does not adversely affect the function of the instrument. To avoid this, all holders and the interior of the autoclave should be checked regularly for corrosion.

4. Release

The preparation must be documented according to regulations and the instruments released for reuse or storage.

5. Storage

Prepared and packaged instruments must be stored at room temperature in a clean, dry place, protected from dust and free of vermin. The period of storage is subject to regulations, the quality of the packaging material, the permeability of sealed seams and storage conditions and generally comprises no more than six months.

Declaration:

The manufacturer has validated that the instructions given above are suitable for preparing a device for reuse. The person performing the preparation is responsible to ensure that the reprocessing actually performed with the equipment used, materials and personnel in the reprocessing facility achieves the desired results. This usually requires validation and routine monitoring of the procedure so that the recommended guidelines from the corresponding institutions (such as the Robert Koch Institute) are complied with. Deviations by the person performing the preparation from the instructions provided should be carefully evaluated for their effectiveness and possible advantageous consequences.

Contact information for the manufacturer

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All rights reserved.
Hopf, Ringleb & Co. GmbH & Cie, Gardeschützenweg 82, 12203 Berlin, Alemania